


UNICKSKIN

YOUR EFFECTIVE SOLUTION

UNICLED NECK & CHIN MASK

MULTI-ACTION TREATMENT WITH LED TECHNOLOGY
LIFTS NECK · REDUCES DOUBLE CHIN · SLIMS JAWLINE

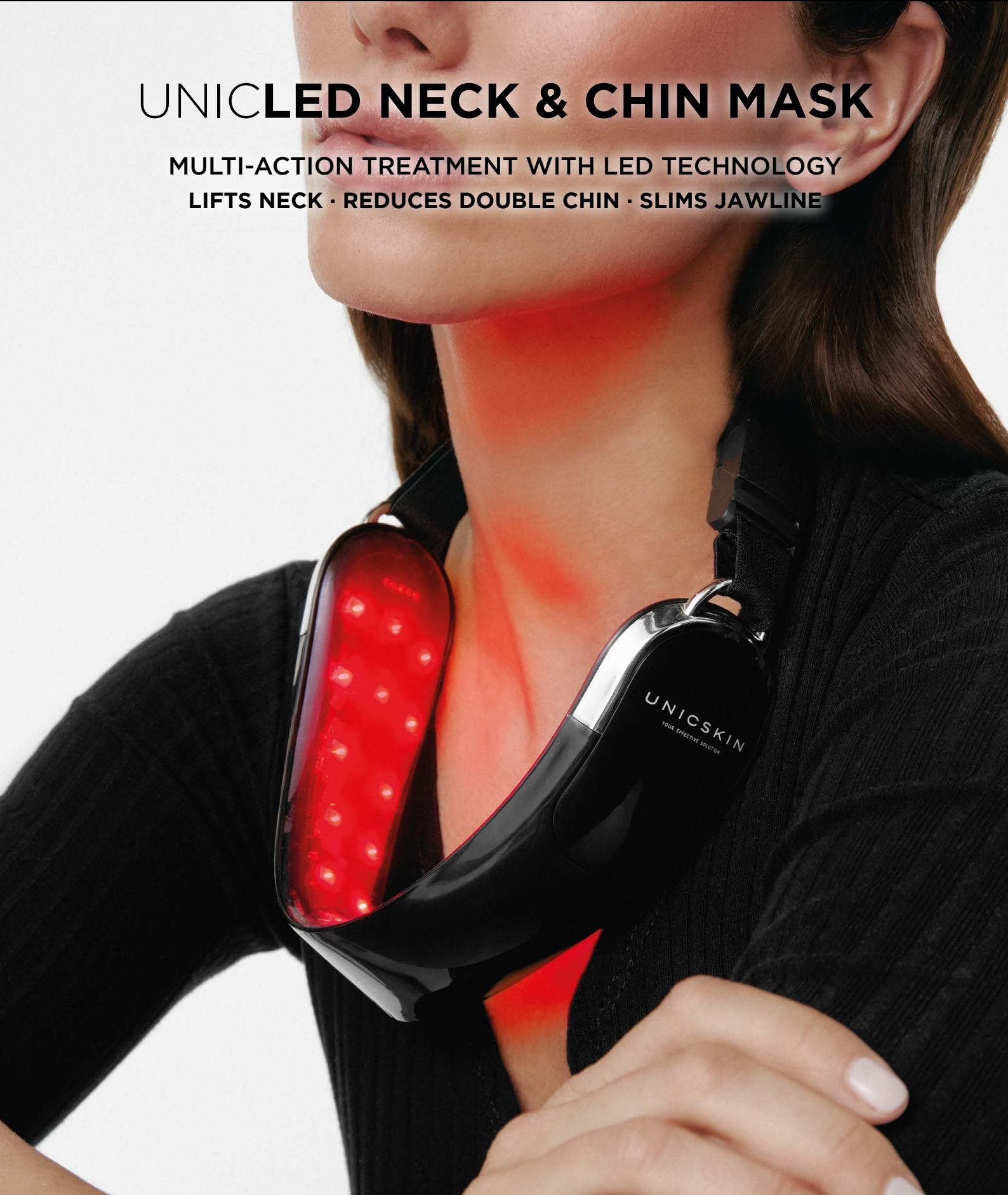

USER MANUAL



LET'S GET STARTED

Congratulations! You have taken the first step toward having a youthful and radiant face, and you have done it with UNICSKIN, global experts in Techno-beauty devices that effectively care for your skin from the comfort of your own home.

Before starting to use this device, please read all instructions and precautions carefully and use it solely for the purpose described in this manual.

LED TECHNOLOGY

UNICLED NECK & CHIN MASK is a diode LED light emitting device that offers a safe combination of different wavelengths. **It acts on the skin and cells, thus reducing oxidative stress and stimulating cellular energy production.** The product offers two modes: neck and chin. Depending on the mode you choose, the device will directly apply one function or another, without you having to do anything. The main benefits that you will notice will be: reduction of neck wrinkles, thinning of the double chin and jawline slimming. In addition, the researchers argue **that LED light therapy attacks adipocytes, which are cells that store fat and cause lipids to disperse. This therapy helps the body eliminate fat cells**, providing a supreme triple action when applied to the neck and chin:


- LIFTS NECK
- REDUCES DOUBLE CHIN
- SLIMS JAWLINE

USE IT 4 TIMES A WEEK. VISIBLE RESULTS IN 4 WEEKS.


WHAT TO EXPECT?

Users typically see a reduction in lines and wrinkles, and a loss of fat within 2 weeks of use (4 times a week). Optimal results will start to be seen after 1 month of treatment.


LED LIGHTS WITH DEEP PENETRATION FOR BRIGHTER AND FIRMER SKIN TEXTURE:

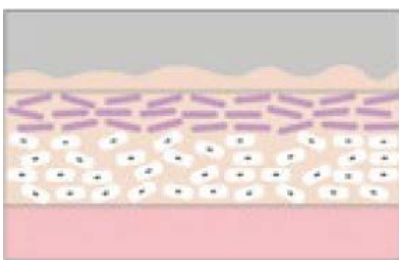

Weaker elasticity.
Deeper wrinkles and fine-lines.

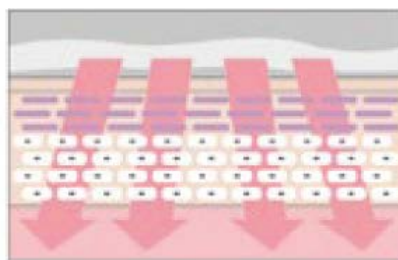

RED light therapy.
Deep penetration into the skin.

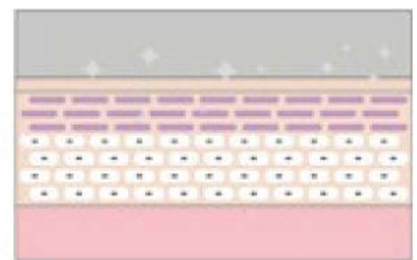

Stronger elasticity.
Smoother wrinkles and fine-lines.



COMPLETE YOUR ROUTINE

We recommend that you combine the mask with the serum that best suits your skin needs.



TIMER

The effectiveness of the treatment includes 20-minute sessions.



CONTROL BUTTONS

Select the mode of use according to the area to be treated: neck, décolleté or double chin.



COMPLEMENTS

Power charging cable.

BENEFITS

**1 DEVICE, 2 MODALITIES, 3 BENEFITS.
EACH MODALITY LASTS 20 MINUTES.**

CHOOSE THE TREATMENT THAT BEST SUITS YOUR NEEDS:

CHIN MODE

DOUBLE CHIN

850nm

Press the power button and choose the "CHIN" mode (you will see a light on above the word). For 20 minutes, the device will operate by:

Infrared LED light (850nm) helps to slim down the double chin and face, defining the jawline or the V-Face.



- ✓ Slims down the face
- ✓ Defines the V-Face
- ✓ Reduces the double chin



Complete your CHIN routine with UNICLIFTING+ firming serum (instant lifting effect).

NECK MODE


NECK & DÉCOLLETÉ

650nm

- ✓ Reduces wrinkles
- ✓ Improves sagging
- ✓ Minimizes signs of aging



Complete your NECK routine with UNICWRINKLE-STOP anti-wrinkle serum (natural botox effect).

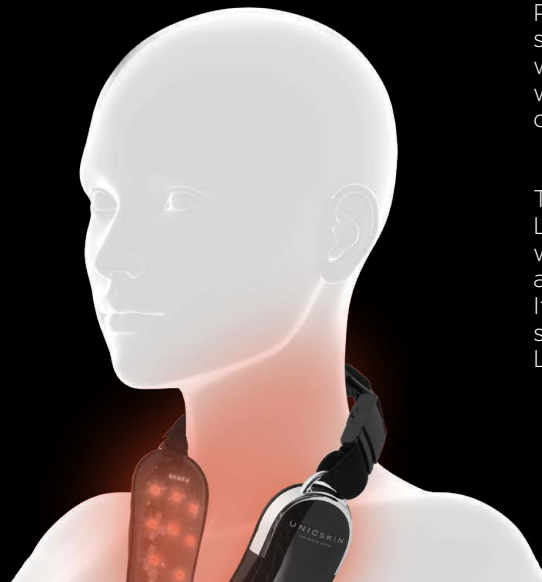

Press the power button and select the "NECK" mode (you will see a light on above the word). For 20 minutes, the device will operate by means of:

The effectiveness of Red LED Light (650nm) that is due to its wavelength, which reaches all areas of the neck and décolleté. Its efficacy covers the entire surface of the skin exposed to LED light emission, as it is shown in the image.

OPERATING INSTRUCTIONS

SET UP


01. Take the UNICLED NECK & CHIN MASK out of the box.
02. Plug the USB charging cable into the device.
03. Plug USB charging cord into outlet for charging.
04. Charge it for three hours (a full charge of the device takes 6 sessions).

HOW TO USE



STEP 1

Before starting the treatment, cleanse the face, neck and décolleté area without applying your daily beauty routine. Turn on the device by pressing and holding the button for 1 second.

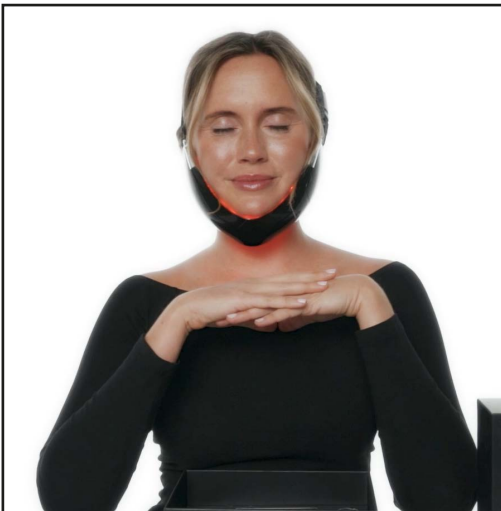


STEP 3

Place the device on the double chin or neck and décolleté and leave it on for 20 minutes. The device will automatically turn off.



STEP 2

Once the mask is turned on, select the mode of use according to the area to be treated: NECK for the neck and décolleté or CHIN for the double chin.

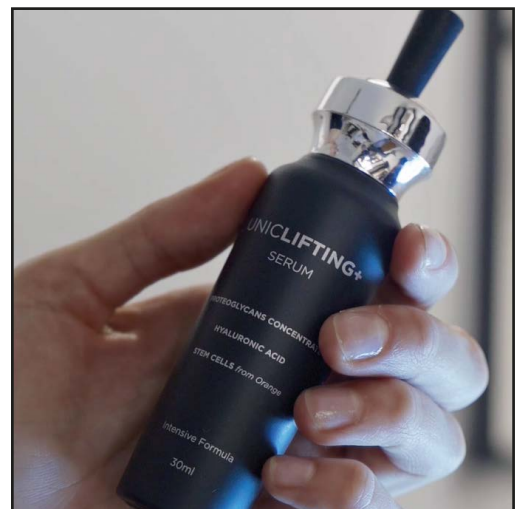

STEP 4

Once the treatment is finished, apply the serum recommended for the mode (see more in the following section: "Benefits"). For optimal results, use the device 3 days a week for 4 weeks. For maintenance, 1-2 times a week.

IMPORTANT SAFETY INFORMATION

- Do not look directly into the lights. The device should be positioned so that it does not cover the eyes.
- Do not leave the device in direct sunlight or in temperatures above 107°F/42°C as this may damage the plastic coating of the device.
- This device is intended for use by adults 18 years of age and older, and is not recommended for use by children under 12 years of age.
- Please consult your physician before using this device if: you are or suspect you are pregnant, nursing, if you have thyroid problems, if you are taking prescription or over-the-counter medications as these cause sensitivity to light.
- Do not place the device near water or any wet surface. If the device comes into contact with water, do not touch it until you unplug it.
- Do not use the device on open wounds.
- Modification of this device is not permitted.
- Always unplug the device from the outlet after charging is complete.
- Do not forcefully bend or abuse the USB charging cable. This will void the warranty.
- Discontinue use of the device immediately if you believe the device has any defects or causes any allergies.

TECHNICAL SPECIFICATIONS

Contents: device (x1), USB-C charger cable (x1) and manual (x1).

Power: 5W

Intensity: Anti-aging 20mW/cm² | Slimming 50mW/cm²

Voltage: 5V

Wavelength: 630nm (red) 660nm (deep red) 850nm (infrared)

Number of lights: 112

DISCLAIMER

This product is not intended to diagnose, cure, or prevent any disease. UNICSKIN S.L. makes no claims, representations, or warranties regarding the ability of this product to cure any mental, physical, or skin condition. A qualified health professional should be consulted with regards to any condition requiring medical attention.

WARRANTY

Except for manufacturing defects, we do not accept returns of the product purchased through our website (www.unicskin.com) once it has been used and put into operation. If the product was purchased outside our website, please refer to the warranty policy offered by the retailer, which in no case shall replace this warranty policy.

If the packaging has not been opened, UNICSKIN will accept the return within 14 days from the date of purchase of the product, at the customer's expense and risk. The device must be returned in the same packaging, together with the original proof of purchase. Please contact Customer Service at the following e-mail address customerservices@unicskin.com, or contact telephone number 912 797 377 to receive information on the address to which the product should be sent.

In case of manufacturing defect and during the first two years from the date of purchase, the return of any of our Techno-Beauty products may result in a replacement for the same product (in which case the cost of return will be borne by UNICSKIN), a refund or a reduction in the price of the product in the event that the defect detected is minor or the value of the good at the time of making the claim does not justify its replacement due to the time elapsed.

It will be necessary in any case to keep and present the receipt or invoice corresponding to the purchased product where the date of purchase is clearly legible.

In addition, UNICSKIN reserves the right to carry out the appropriate checks on the product purchased, requesting photographs or evidence of the alleged non-conformity or the shipment of the defective product, with the shipping costs initially being borne by the customer, and proceeding to reimbursement once it has been verified by UNICSKIN that the fault detected corresponds to a manufacturing defect.

Excluded from this warranty policy are those products that have suffered damage due to storage conditions, maintenance, or any other type of handling by the customer. The product complies with the legal and regulatory standards applicable in the European Community at the date of manufacture and distribution.

CLINICAL RESEARCHS

Conducted by the U.S. Department of Health and Human Services (U.S. National Library of Medicine).

In a study published in Lasers in Surgery Medicine, participants who received light therapy showed a significant reduction in overall circumference, including numerous parts of the body. The researchers concluded that light therapy can reduce the circumference of specific areas of the body that are treated with natural light. (1)

A study conducted by Brazilian researchers in 2015 evaluated the effects of light therapy and exercise in 64 obese women, aged 20 to 40 years. One group exercised and received red light therapy treatments, while the control group only exercised. They concluded that the combination of red light therapy and exercise were more effective than exercise alone in reducing fat mass (2).

Red light therapy is safe, non-invasive, uses no chemicals or drugs, and is not associated with harmful side effects. Red light therapy simply harnesses the natural healing and rejuvenating benefits of a specific range of therapeutic natural light and delivers this specific energy at a higher rate than the sun, without the damage of UV rays. A major 2013 study on light therapy by Harvard and MIT researchers praised light therapy for its "non-invasive nature and its near-total absence of side effects." (3)

(1) Jackson RF, Dedo DD, Roche GC, et al. Low-level laser therapy as a noninvasive approach to body contouring: a randomized, controlled study. Lasers in Surgery and Medicine. December 2009; 41 (10): 99-809. www.ncbi.nlm.nih.gov/pubmed/20014253.

(2) Sene-Fiores M, Duarte FO, et al. The potential of phototherapy to reduce cor-poral fat, insulin resistance, and obesity-related "metabolic inflexibility" in women undergoing weight loss treatment. Laser Surg Med. October 2015; 47 (8): 634-42. <https://www.ncbi.nlm.nih.gov/pubmed/26220050>.

(3) Avci P, Gupta A, Sadasivam M et al. "Low level laser therapy (LLLT) on the skin: stimulating, healing, restorative". Semin Cutan Med Surg. March 2013; 32 (1): 41-52. <https://www.ncbi.nlm.nih.gov/pmc/articles/PMC4126803/>.

ACCESSORIES UNICLED NECK & CHIN MASK



User
manual



USB charger



Bag

CUSTOMER SERVICE

If you have any questions or comments you would like to share with us, please write to us at the following email address: customerservices@unicskin.com.

UNICSKIN

YOUR EFFECTIVE SOLUTION

GLOBAL EXPERTS IN TECHNOBEAUTY & BEAUTY AT HOME
EFFECTIVE RESULTS PROVEN BY SCIENCE
SUSTAINABILITY & FORMULAS +94% OF NATURAL ORIGIN
INNOVATIVE LUXURY AT YOUR FINGERTIPS
